



Code
2979




Total Chlorine ppm



Using a cup-size sample, immerse strip & **SWIRL VIGOROUSLY 3 TIMES**. Remove strip with **pad face up**. **DO NOT SHAKE OFF EXCESS WATER**. With pad still facing up, wait **10 seconds** and **read** the strip. Discard strip after 30 seconds.

TIPS

- Keep wet fingers out of the bottle. Cap bottle tightly.
- Store in a cool, dry place.
- For large vats, swirl one time only.

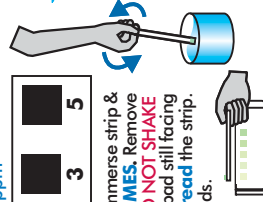

Lot # 7126 Exp. 8/09

Chelsterow, MD 21620 www.lamotte.com/watcon